



Sanyo HEC-DR6700K - Zero Gravity Massage Chair Instruction Manual

Massage lounger

1
2
3
4
5
6
7
8
9
10
11
12
13
14
15
16
Table Of Contents
17



-

[Table of Contents](#)

-

Troubleshooting

•

Bookmarks

Quick Links


[1 Remote Control](#)

[2 Maintenance Procedures](#)

[3 Troubleshooting](#)

[4 Specifications](#)

[Download this manual](#)

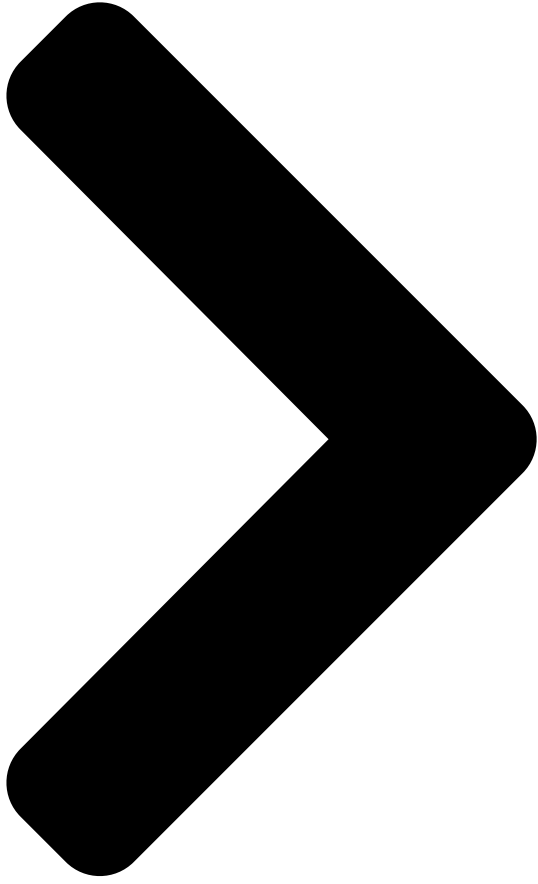

SANYO MESSAGE LOUNGER LIMITED WARRANTY

OBLIGATIONS

In order to obtain warranty service, the product must be delivered to and picked up from an Authorized Sanyo Factory Service Center at the user's expense, unless specifically stated otherwise in this warranty.

If you need assistance with product operation, authorized service center referral, service assistance, problem resolution, or parts and accessories availability; you may contact us by:

A. Visiting www.sanyoservice.com on a 24 hr/7 days a week basis or
 B. calling the customer support center at 1-800-421-5013
 Weekdays 7:30 AM - 5:00 PM Pacific Time

SANYO MESSAGE LOUNGER LIMITED WARRANTY

(USA ONLY)

OBLIGATIONS

THIS WARRANTY IS VALID ONLY ON SANYO PRODUCTS PURCHASED AND USED IN THE UNITED STATES OF AMERICA. THIS WARRANTY APPLIES ONLY TO THE ORIGINAL RETAIL USER AND DOES NOT APPLY TO PRODUCTS USED FOR ANY INDUSTRIAL, PROFESSIONAL OR COMMERCIAL PURPOSE. THE ORIGINAL DATED BILL OF SALE OR SALES SLIP MUST BE SUBMITTED TO THE AUTHORIZED SANYO SERVICE CENTER AT THE TIME WARRANTY SERVICE IS REQUESTED.

Subject to the OBLIGATIONS above and EXCLUSIONS below, SANYO FISHER COMPANY (SFC) warrants this SANYO product against defects in materials and workmanship for the periods of LABOR and PARTS specified below. SFC will repair or replace (at its option) the product and any of its parts which fail to conform to this warranty.

The warranty period commences on the date the product was first purchased at retail.

LABOR	PARTS	*ROLLER MECHANISM	**IN - HOME
1 YEAR	1 YEAR	3 YEARS	1 YEAR

Subject to the OBLIGATIONS above and EXCLUSIONS below, SANYO FISHER COMPANY (SFC) warrants this SANYO product against defects in materials and workmanship for the periods of LABOR and PARTS specified below. SFC will repair or replace (at its option) the product and any of its parts which fail to conform to this warranty.

The warranty period commences on the date the product was first purchased at retail.

Subject to the OBLIGATIONS above and EXCLUSIONS below, SANYO FISHER COMPANY (SFC) warrants this SANYO product against defects in materials and workmanship for the periods of LABOR and PARTS specified below. SFC will repair or replace (at its option) the product and any of its parts which fail to conform to this warranty.

The warranty period commences on the date the product was first purchased at retail.

Subject to the OBLIGATIONS above and EXCLUSIONS below, SANYO FISHER COMPANY (SFC) warrants this SANYO product against defects in materials and workmanship for the periods of LABOR and PARTS specified below. SFC will repair or replace (at its option) the product and any of its parts which fail to conform to this warranty.

The warranty period commences on the date the product was first purchased at retail.

Note: The roller mechanism consists of the mechanical parts directly behind the back rest of the massage chair, excluding belts.

An in-home service is defined as product service in the consumer home or offsite if further repairs are required.

LABOR	PARTS	*ROLLER MECHANISM	**IN - HOME
1 YEAR	1 YEAR	3 YEARS	1 YEAR

This warranty does not cover (A) the adjustment of customer-operated controls as explained in the appropriate model's instruction manual, or (B) the repair of any problem whose solution has been identified in the instruction manual.

This warranty shall not apply to the following items: non-electrical/mechanical attachments; accessories; damage to head cover, back pad, arm covers, seat cover, foot rest cover when not caused by defects in workmanship.

This warranty shall not apply to setup, installation, removal of the product for repair, or reinstallation of the product after repair. This warranty does not apply to repairs or replacements necessitated by any cause beyond the control of SFC including, but not limited to, any malfunction, defect or failure caused by or resulting from unauthorized service or parts, improper maintenance, operation contrary to furnished instructions, shipping or transit accidents, modification or repair by the user, abuse, misuse, neglect, accident, incorrect power line voltage, fire, flood or other Acts of God, or normal wear and tear.

EXCLUSIONS

The foregoing is in lieu of all other expressed warranties and SFC does not assume or authorize any party to assume for it any other obligation or liability.

THE DURATION OF ANY WARRANTIES WHICH MAY BE IMPLIED BY LAW (INCLUDING THE WARRANTIES OF MERCHANTABILITY AND FITNESS) IS LIMITED TO THE TERM OF THIS WARRANTY. IN NO EVENT SHALL SFC BE LIABLE FOR SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING FROM OWNERSHIP OR USE OF THIS PRODUCT, OR FOR ANY DELAY IN THE PERFORMANCE OF ITS OBLIGATIONS UNDER THIS WARRANTY DUE TO CAUSES BEYOND ITS CONTROL.

SOME STATES DO NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS AND/OR DO NOT ALLOW THE EXCLUSION OR LIMITATION OF CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATIONS AND EXCLUSIONS MAY NOT APPLY TO YOU.

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS. YOU MAY HAVE OTHER RIGHTS, WHICH VARY FROM STATE TO STATE. For your protection in the event of theft or loss of this product, please fill in the information below for your own personal records.

Model No. _____ Serial No. _____
 (Located on back or bottom side of unit)

Date of Purchase _____ Purchase Price _____

Where Purchased _____

EXCLUSIONS

Serial No. _____
 (Located on back or bottom side of unit)

Purchase Price _____

Please register your message chair at www.sanyomassagechairs.com

98F-BP111336CB-0 0308

(USA ONLY)
 An ROLLER MECHANISM
 An IN - HOME
 3 YEARS
 1 YEAR

Where Purchased _____

Serial No. _____
 (Located on back or bottom side of unit)

Purchase Price _____



SANYO

SANYO Fisher Company
 21605 Plummer Street
 Chatsworth, CA 91311
www.sanyo.com

**INSTRUCTION MA
 Massage Lounger**



Be sure to read this manual thoroughly in order to ensure proper use. In particular, make "IMPORTANT SAFETY INSTRUCTIONS". Keep this manual in a safe place so that it can be found easily.

**INSTRUCTION MANUAL
 Massage Lounger**

Be sure to read this manual thoroughly in order to ensure proper use. In particular, make certain to read the section "IMPORTANT SAFETY INSTRUCTIONS".

Keep this manual in a safe place so that it can be found easily.

SANYO MESSAGE LOUNGER LIMITED WARRANTY

(USA ONLY)

OBLIGATIONS

THIS WARRANTY IS VALID ONLY ON SANYO PRODUCTS PURCHASED AND USED IN THE UNITED STATES OF AMERICA. THIS WARRANTY APPLIES ONLY TO THE ORIGINAL RETAIL USER AND DOES NOT APPLY TO PRODUCTS USED FOR ANY INDUSTRIAL, PROFESSIONAL OR COMMERCIAL PURPOSE. THE ORIGINAL DATED BILL OF SALE OR SALES SLIP MUST BE SUBMITTED TO THE AUTHORIZED SANYO SERVICE CENTER AT THE TIME WARRANTY SERVICE IS REQUESTED.

Subject to the OBLIGATIONS above and EXCLUSIONS below, SANYO FISHER COMPANY (SFC) warrants this SANYO product against defects in materials and workmanship for the periods of LABOR and PARTS specified below. SFC will repair or replace (at its option) the product and any of its parts which fail to conform to this warranty.

The warranty period commences on the date the product was first purchased at retail.

LABOR	PARTS	*ROLLER MECHANISM	**IN - HOME
1 YEAR	1 YEAR	3 YEARS	1 YEAR

Subject to the OBLIGATIONS above and EXCLUSIONS below, SANYO FISHER COMPANY (SFC) warrants this SANYO product against defects in materials and workmanship for the periods of LABOR and PARTS specified below. SFC will repair or replace (at its option) the product and any of its parts which fail to conform to this warranty.

The warranty period commences on the date the product was first purchased at retail.

If you need assistance with product operation, authorized service center referral, service assistance, problem resolution, or parts and accessories availability; you may contact us by:

A. Visiting www.sanyoservice.com on a 24 hr/7 days a week basis or
 B. calling the customer support center at 1-800-421-5013
 Weekdays 7:30 AM - 5:00 PM Pacific Time

CONTENTS

PAGE

IMPORTANT SAFETY INSTRUCTIONS Ac Ac 1-4

NAME/FUNCTION OF EACH PART Ac Ac 5-8

HOW TO PLACE THE MASSAGER Ac Ac 9

HOW TO INSTALL THE STAND Ac Ac 10-11

HOW TO USE Ac Ac 12-13

PREPARATIONS PRIOR TO USE Ac Ac 14-15

ADJUSTMENT OF CHAIR Ac Ac 16-18

MESSAGING Ac Ac 19-21

WHOLE BODY SENSOR AUTOMATIC COURSES Ac Ac 16-18

AUTOMATIC COURSES Ac Ac Ac Ac Ac Ac 19-21

MANUAL COURSES Ac Ac Ac Ac Ac Ac Ac Ac 21-24

ADJUSTING MESSAGE INTENSITY Ac Ac Ac 25

AFTER USE Ac Ac Ac Ac Ac Ac Ac Ac Ac Ac Ac 26

MAINTENANCE PROCEDURES Ac Ac Ac Ac Ac 26

CONTENTS OF MASSAGING Ac Ac Ac Ac Ac 27

TROUBLE-SHOOTING Ac Ac Ac Ac Ac Ac Ac Ac 28

SPECIFICATIONS Ac Ac Ac Ac Ac Ac Ac Ac 29

**INSTRUCTION MA
 Massage Lounger**

[Next Page](#)

1
2
3
4
5

Related Manuals for Sanyo HEC-DR6700K - Zero Gravity Massage Chair

[Massager Sanyo HEC-DR6700K - Zero Gravity Massage Chair Parts List](#)

Massage lounger (17 pages)

[Massager Sanyo HEC-DR7700C Instruction Manual](#)

Massage lounger (16 pages)

[Massager SANYO HEC-DR7700 Service Manual](#)

Chair type massager(massage lounger) (29 pages)

[Massager Sanyo HEC-DR7700 Service Manual](#)

(28 pages)

[Massager SANYO HEC-DR7700 Instruction Manual](#)

Chair type massager (17 pages)

[Massager SANYO HEC-DR8700 Instruction Manual](#)

Massage lounger (36 pages)

[Massager Sanyo HEC-DR8700 Service Manual](#)

Chair type massager (massage lounger) (35 pages)

[Massager SANYO HEC-DR5000 Service Manual](#)

Chair type massager (30 pages)

[Massager Sanyo HECSA5000 Service Manual](#)

Chair type massager (massager lounger) (31 pages)

[Massager Sanyo HECSR1000K - Stiffness Sensor - Multi Roller Massage Chair Service Manual](#)

Service manual (20 pages)

[Massager Sanyo HEC-DR5700 Instruction Manual](#)

(32 pages)

[Massager Sanyo HEC-SR1000K Instruction Manual](#)

Sanyo hec-sr1000k massager: user guide (24 pages)

[Massager Sanyo HEC-RX1BK Instruction Manual](#)

Relaxation chair (9 pages)

[Massager Sanyo HEC-SA5000K Instruction Manual](#)

Sanyo hec-sa5000k massager: user guide (28 pages)

[Massager SANYO HEC-A2000 Service Manual](#)

Chair type massager (19 pages)

[Massager Sanyo HEC-A3700 Instruction Manual](#)

Massage lounger (24 pages)

Summary of Contents for Sanyo HEC-DR6700K - Zero Gravity Massage Chair

[Page 1](#) SANYO MESSAGE LOUNGER LIMITED WARRANTY (USA ONLY) OBLIGATIONS In order to obtain warranty service, the product must be delivered to and picked up from an Authorized Sanyo Factory Service INSTRUCTION MANUAL Center at the user's expense, unless specifically stated otherwise in this warranty.

[Page 2: Important Safety Instructions](#)

IMPORTANT SAFETY INSTRUCTIONS All of the following information is important. Be sure to follow the instructions carefully. When using an electrical appliance, basic precautions should always be followed, including the Do not wind the power cord, remote control cord or sensor controller cord around the main unit, following: remote control, sensor controller, or remote control stand.

Page 3 **IMPORTANT SAFETY INSTRUCTIONS WARNING** Do not attempt to disconnect the power plug by pulling the cord, but always by holding the power plug and pulling it out. Also, clean the dust off the power plug periodically. Otherwise, electric shock, short circuit or fire may occur. Be careful not to spill water over the remote control, sensor controller, seat, or foot rest.

[Page 4: Name/Function Of Each Part](#)

NAME/FUNCTION OF EACH PART Remote control HOME POSITION button Used to return the back rest, foot rest and seat to Head cover their home positions. All display icons are shown here Seat for explanation. ZERO-G button Used to adjust the back rest, foot rest and seat to Sensor controller the preset zero gravity positions.

[Page 5: How To Place The Massager](#)

NAME/FUNCTION OF EACH PART HOW TO PLACE THE MASSAGER Positions that can be massaged. PRESSURE display How to move the massager Positions where stiffness has been detected. Displays the strength of the pressure placed on the massaging rollers (the pressure from the body). Raise the back rest.

[Page 6: How To Install The Stand](#)

HOW TO INSTALL THE STAND Installation on the right-hand side of the chair Installation on the left-hand side of the chair Tools needed: Phillips (crosshead) screwdriver, flathead screwdriver (All screwdrivers are not included) Tools needed: Phillips (crosshead) screwdriver, flathead screwdriver (All screwdrivers are not included) Remove the cover of the right-hand Put the remote control on the Remove the covers on the outer...

[Page 7: How To Use](#)

HOW TO USE PREPARATIONS PRIOR TO USE The procedure for use of the massager is explained here. See the referred pages for details. **WARNING** Check the cloth of the chair. Prior to use, be sure to raise the back pad The following persons are requested to consult a Prior to use, make certain to read the "IMPORTANT SAFETY INSTRUCTIONS".

[Page 8: Adjustment Of Chair](#)

ADJUSTMENT OF CHAIR Before sitting on the massage chair Press ZERO-G button to move the back rest, foot rest, and the seat to achieve the desired Bring all parts of the massage chair to reclining position. Foot rest Press their home position by pressing the HOME POSITION button.

Page 9 **ADJUSTMENT OF CHAIR MASSAGING (WHOLE BODY SENSOR AUTOMATIC COURSES)** Adjust the angle of the foot rest. Press the OPERATION ON/OFF button. **WARNING** When the button is pressed, the lamp of the remote control Foot rest can be adjusted independently. illuminates and the sole warmer will turn on. When the foot rest is suddenly released Keep the FOOT REST button pressed until When the foot rest is fully lowered, the foot rest will rise a little.

Page 10 **MASSAGING (WHOLE BODY SENSOR AUTOMATIC COURSES)** Measurement of shoulder position Massaging operation and display of remote control Starts massaging automatically. Displays the aim of massaging positions, the type Measure the position of shoulder and Measure the position of shoulder. of massaging, the stiffness, body pressure and the match the position of massaging rollers Check if the position of shoulder is correct.

Page 11 **MASSAGING (AUTOMATIC COURSES)** Measurement of shoulder position Massaging operation and display of remote control Press the OPERATION ON/OFF button. Displays the aim of massaging positions, the type When the button is pressed, the lamp of the remote Measure the position of shoulder and of massaging, body pressure and the remaining control illuminates and the sole warmer will turn on.

Page 12 **MASSAGING (MANUAL COURSES)** Massaging of upper body Press the MANUAL button. MANUAL button Press the desired MANUAL (UPPER BODY) MESSAGE button. Refer to page 23 for details of the contents of massage. For grasping shoulder massaging and spine

stretching, shoulder position is first measured.

[Page 13](#) MASSAGING (MANUAL COURSES) Kinds of massages and adjustments of massaging rollers
Massaging of lower body Several types of desired massage can be selected.
Simultaneous massaging of any part of "upper body" and "lower body (air)" is possible. It is possible to adjust the movement speed, height position and width of massaging rollers and the strength of massage. Independent massaging of "upper body"...

[Page 14: Adjusting Massage Intensity](#)

ADJUSTING MESSAGE INTENSITY AFTER USE CAUTION When the massage is too strong After use, be sure to set the POWER switch to the [OFF] position, and the LOCK switch to the Use the massager after placing [LOCK] position. When the equipment will not be used for an extended period, be sure to disconnect the power a towel, etc., on the back pad.

[Page 15: Contents Of Massaging](#)

∅ñ The drawing of the massage for each course is omitted in part, in order to describe the CONTENTS OF MASSAGING operation in a simplified manner. Also, the drawings of air massage operation for lower body are all omitted. WHOLE BODY SENSOR AUTOMATIC COURSE WHOLE BODY SENSOR AUTOMATIC COURSE ∅q ∅r...

[Page 16: Troubleshooting](#)

TROUBLE-SHOOTING SPECIFICATIONS WARNING Power consumption 290 W If the trouble still persists even after the following proper countermeasures have been taken, stop using the massager to avoid possible accidents. Disconnect the power plug from the outlet and consult an authorized service center. Power dissipation of electric heating equipment 33 W Do not attempt to repair the equipment at home by yourself, as it is very dangerous.